



DISPOSABLE SYRINGE

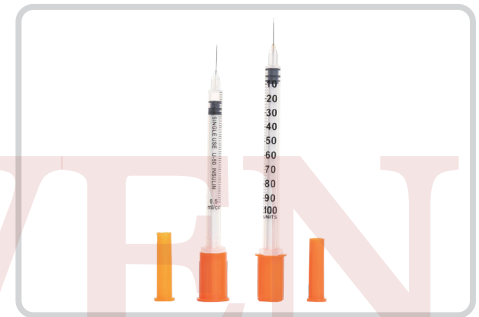


Feature:

- Sterilized by EO gas.non-toxic. non-pyrogenic,single use only.
- Luer lock or luer slip.
- Sufficiently transparent barrel allows easy measurement.
- The plunger fits the inside of barrel very well to allow for free and smooth movement.
- Material for gasket: Natural latex/ Synthetic rubber (Latex free).
- Material for barrel and plunger: Medical grade PP.
- With or without needle.
- Individual package: blister bag/ PE bag.
- Available size: see the table bellow.

SIZE	SPECIFICATION
1ml	LUER SLIP or LUER LOCK
2ml	LUER SLIP or LUER LOCK
3ml	LUER SLIP or LUER LOCK
5ml	LUER SLIP or LUER LOCK
10ml	LUER SLIP or LUER LOCK
20ml	LUER SLIP or LUER LOCK
30ml	LUER SLIP or LUER LOCK
50ml	LUER SLIP or LUER LOCK
60ml	LUER SLIP or LUER LOCK
60ml	CATHETER

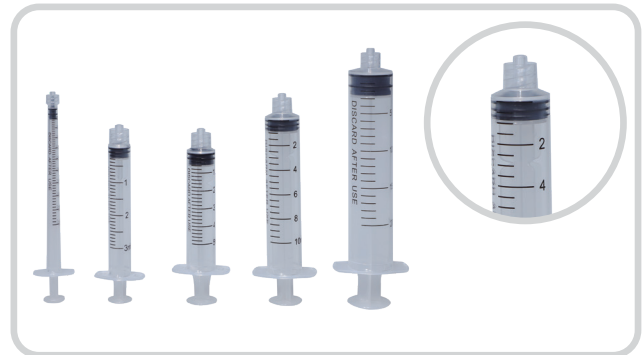


SAFETY RETRACTABLE SYRINGE

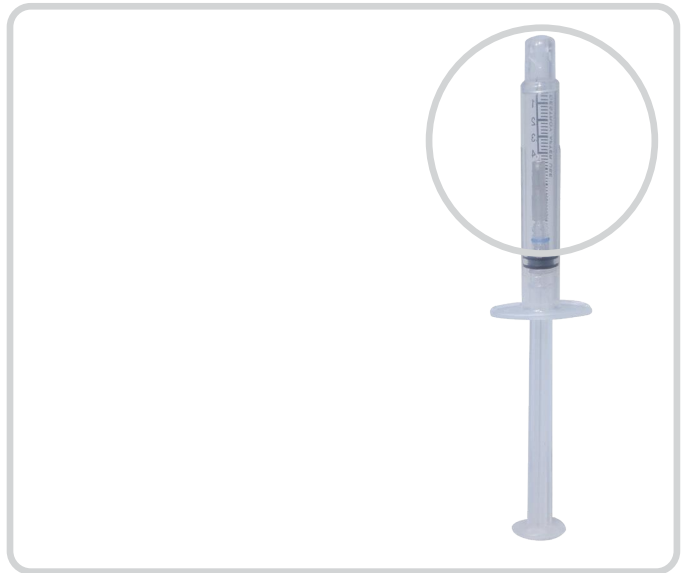

Feature:


- Innovative mechanical secure-locking system design;
- Manual pull-back, simple operation as disposable syringe; Completely eliminate the second use of syringes and the risk of stabbing medical staff;
- The needle retracts into the barrel, safe; No spare parts, avoid unnecessary material pollution and waste;
- No gap;
- National patent.

Using instruction:

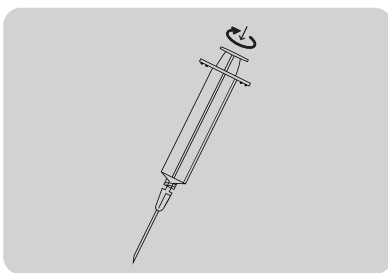


Step1

Push plunger firmly against base of barrel and turn 30° clockwise

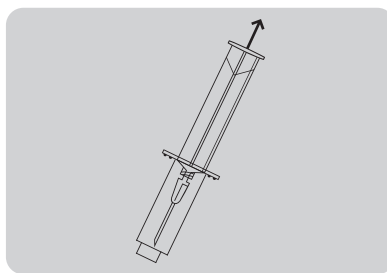

Step2



Pull back plunger to retract needle

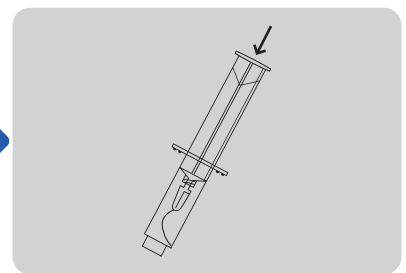

Step3

Push plunger forward to destroy needle


AUTO-DESTROY SYRINGE

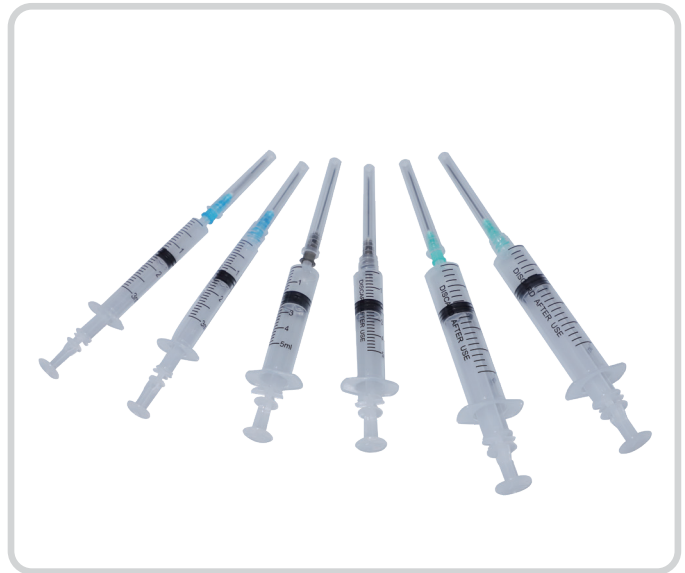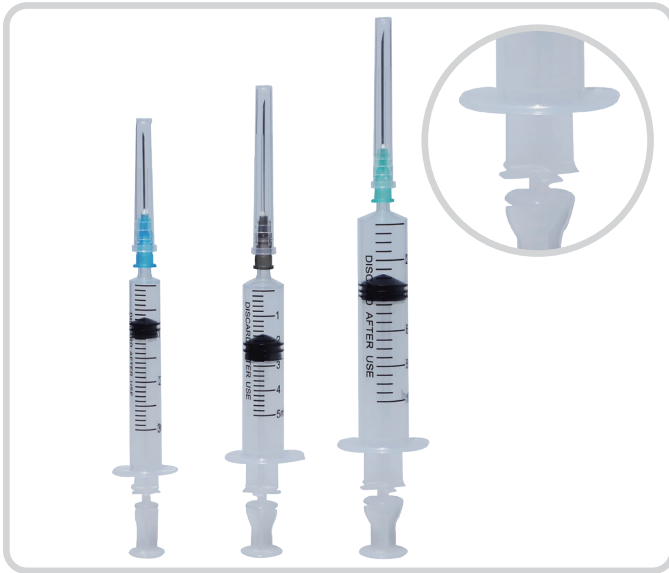

Feature:

- On the basis of traditional syringe, it has functions of auto-destroy, special structure, and easy operation;
- Can be destroyed automatically after single use;
- Easy and safe characteristics, which makes it a good choice for liquid medicine injection;
- Size: 1ml, 2ml, 3ml, 5ml, 10ml;
- Different packages : according to customer requirements.

Using instruction:

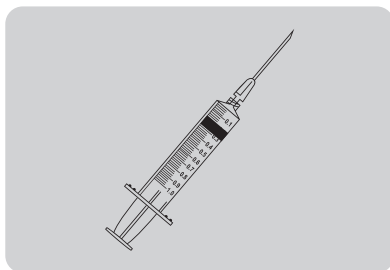


Step1

Pull back the plunger to charge the syringe with solution.



Step2

Press plunger forward to complete injection until it reaches the stop position. Lock mechanism will be activated and locks plunger in stop position. An audible "Click" will sound.

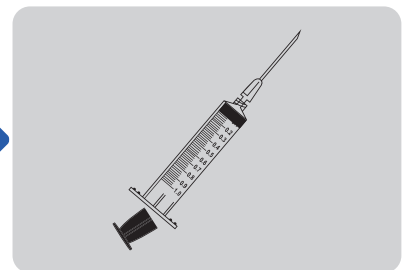

Step3

Forcing the plunger backward will cause it to break, Dispose safely in a disposable container.


INFUSION SET




- Plastic spike with air vent.
- PVC tubing.
- Drip chamber with filter.
- Flow regulator.
- Luer lock connector.



- Plastic spike with air vent.
- PVC tubing.
- Drip chamber with filter.
- Flow regulator.
- Luer lock connector with hypodermic needle.
- Bulb (latex free or latex).

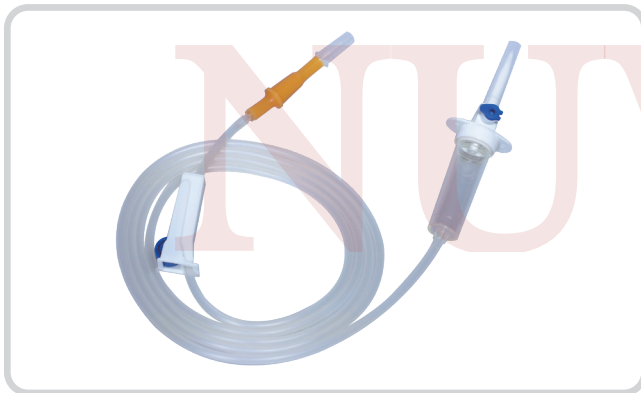


- Plastic spike with air vent.
- PVC tubing.
- Drip chamber without filter.
- Flow regulator.
- Bulb (latex or latex free).
- Luer Slip connector.

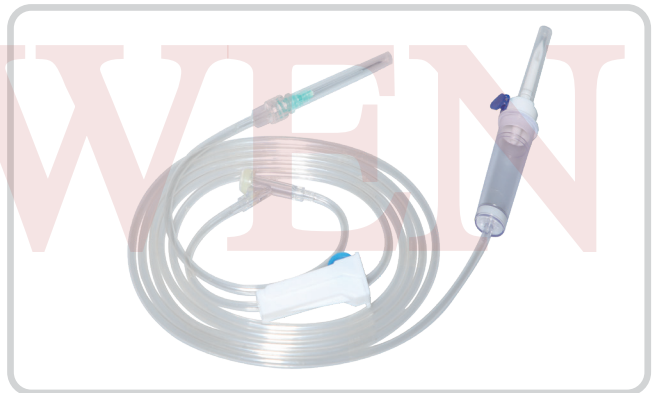

- Plastic spike with air vent.
- PVC tubing.
- Drip chamber with filter.
- Flow regulator.
- Y connector.
- Needle.



- Plastic spike with air vent.
- PVC tubing.
- Drip chamber without filter.
- Flow regulator.
- Y type connector.
- Luer Slip connector.

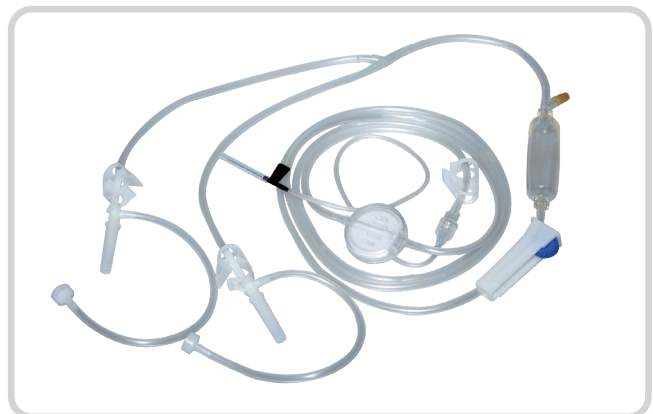



- Drip chamber.
- Flow regulator.
- PVC tubing.
- Liquid filter.
- Robert clamp.
- Scalp vein set.

INFUSION SET




We can customize various types of infusion sets to meet the needs of different customers.


- Small drip chamber 14.3 x 50mm
Medium drip chamber 15.5 x 54mm
Big drip chamber 17.2 x 61mm
- With air vent, without air vent
- Transfer set
- 15drops/ml, 20drops/ml, 60drops/ml
- With filter or without filter
- Small flow regulator, Medium flow regulator, ABS Adelberg flow regulator
- Accurate flow regulator rate of 5ml/hr to a maximum of 300ml/hr
- Latex straight tube, injection bulb, Y injection port.
- Latex free straight tube, injection bulb, Y injection port.
- Needleless Y injection port, Needle free injection cap.
- Luer slip, luer lock connector
- Red cap, blue cap, white cap
- PVC tube length 1.3m to 2m or any required length
- PVC tube DEHP FREE
- Non-PVC tube
- Customization mold construction for any kind of components

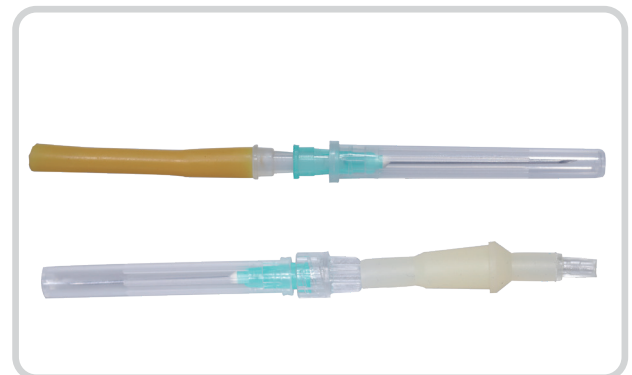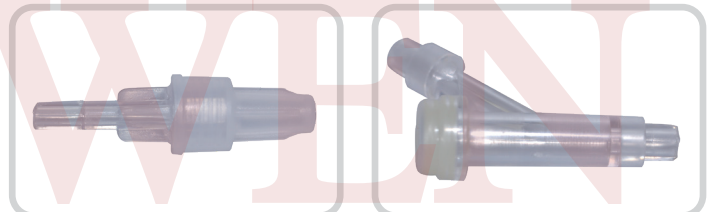


VAGINAL SPECULUMS

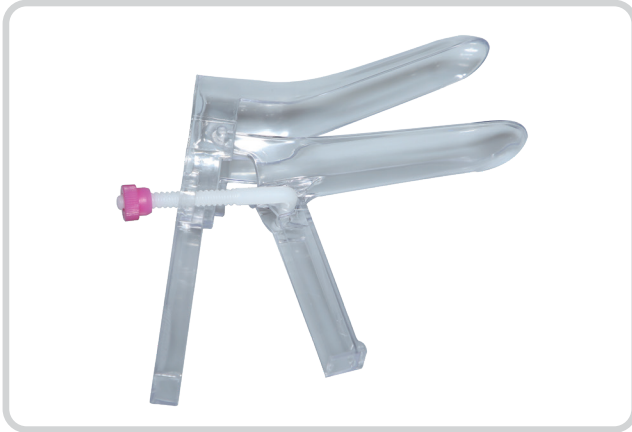

• Side Screw Type
• SIZE:S/M/L

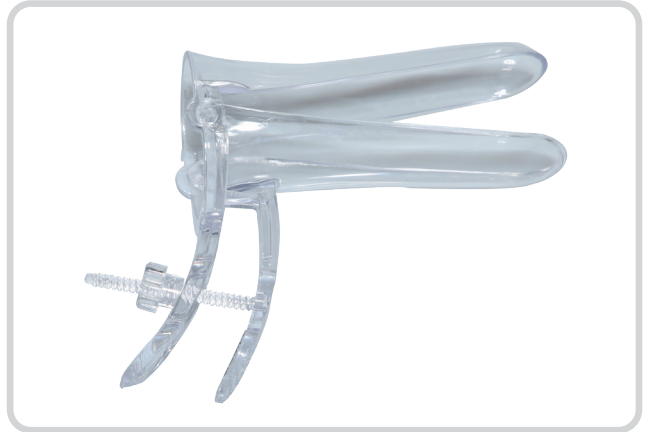

• Middle Screw Type
• SIZE:XS/S/M/L

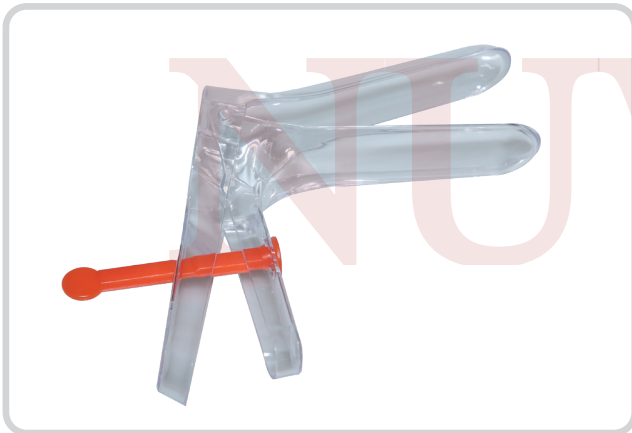

• French Type 1
• SIZE:XXS/XS/S/M/L/XL



• French Type 2
• SIZE:S/M/L

Feature:

- Upper and lower part made of polystyrene.
- Bolt made of polyethylene.
- Perfectly smooth and streamlined edges.
- Color coded bolts.
- For use no more than 60 minutes.
- Single-use.
- Sterile : E.O. Gas.
- Package : PE bag/ Complex film/ Paper-plastic bag/ Bulk.



• Spanish Type
• SIZE:S/M/L

INFUSION SET WITH BURETTE



100ml

150ml

- Air-filter**
- Soft tube**
- Drip chamber**
- Flow regulator**
- Fluid filter**
- Float valve**
- Y injection site (latex or latex free)**

Feature:


- The burette are highly transparent, clear and accurate in graduation, and this makes more convenient operation when infusion.
- Flexible, kink resistant and Non-distortion tubes.
- Credible and agile flow regulator for better flow rate control.
- Vented with hydrophobic bacteria filtration membrane.

Specification:

- Application: Used for the children.
- Type: 100ml, 150ml bottle size, 60 drops/ml.
- Tube length: 150cm.
- Product Instruction: The products are made by medical polymer materials.
- Sterilized by EO gas, non-toxic, pyrogen free.

+86-18680577977

GYNECOLOGICAL PRODUCTS

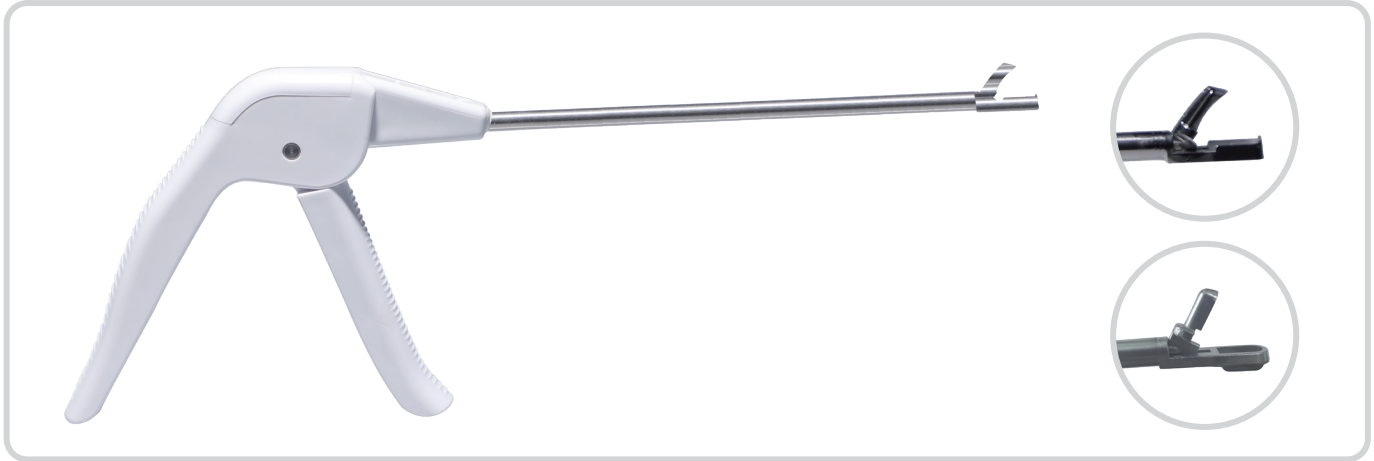

Disposable uterus biopsy forceps

Specification:

- Single use only, sterilized by EO.
- Head material : Stainless steel.
- More Safe, effective and reliable.
- Completely eliminate cross infection in hospital, effectively reduce bleeding of biopsy and avoid pain.
- More medical value.
- Packing: 1pcs/hard blister, 10pcs/box, 60pcs/carton.

















Gynecological set

Specification:

- Latex gloves, L.
- Vaginal Speculums.
- Non-woven fabric, 45 x 45cm, 40 x 60cm.
- The cervical smear.
- The cervical brush.
- Glass slide.

HYPODERMIC NEEDLES

SPECIFICATION				
Size		Color		Length Scope
OD(mm)	Gauge	Color&Code		
0.30	30G		Light yellow	3/8"-11/2"(10mm-40mm)
0.33	29G		Red	
0.36	28G		Blue-green	
0.40	27G		Grey	
0.45	26G		Brown	
0.50	25G		Orange	
0.55	24G		Purple	
0.60	23G		Blue	
0.70	22G		Black	
0.80	21G		Green	
0.90	20G		Yellow	
1.10	19G		Beige	
1.20	18G		Pink	
1.60	16G		White	

Feature:

- Made of high quality stainless steel.
- Semi-transparent needle-hub for the convenient observation of blood reflux.
- Delicately-designed needle-point with features of sharpness, smooth puncture, less damage to the tissue, less painful feeling to the patient.
- Size distinguished by the needle-hub's color for the clear recognition.
- The flexible supply form: in bulk or in blister, sterilized or unsterilized.



Butterfly scalp vein set

Luer Slip and Luer Lock

+86-18680577977

FACE MASK



MEDICAL FACE MASK


Specification:

- Materials : 70% spunbond fabric, 30% melt-blown fabric.
- Standard : EN14683-2019.
- Type : Type I, Type II, Type IIR.
- Size: Adult, Children (17.5 x 9.5cm, 14.5 x 9cm).
- Colour : White, Black, Blue, Floral.

FFP2 MASK

Specification:

- Materials : Non-woven fabric 40%, Melt-blown fabric 30%, Hot air cotton 30%
- Standard : EN149:2001 + A1: 2009 FFP2 NR.
- Size: Adult (16.5 x 10.6cm).
- Colour : White.

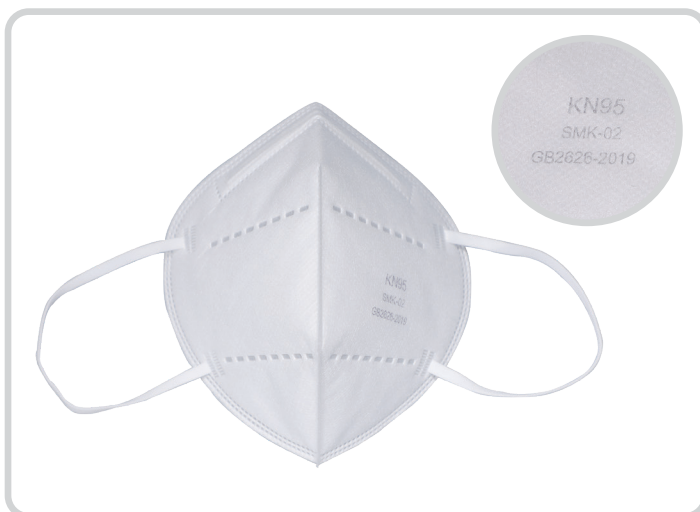

KN95 MASK

Specification:

- Materials : Non-woven fabric 40%, Melt-blown fabric 30%, Hot air cotton 30%
- Standard : GB2626-2019
- Size: Adult (16.5 x 10.6cm)
- Colour : White

URINE BAG



SUPER URINE BAG



SUPER URINE BAG

Feature:

- Made from medical grade PVC .
- 2000ml, with anti-flow device, cross valve.
- With injection port, sampling port.
- Usage : used for collecting bed patient' s urine.



**URINE BAG
SINGLE VALVE**





**URINE BAG
PULL PUSH VALVE**



**URINE BAG
CROSS VALVE**

DISPOSABLE IN VIVO SAMPLING FORCEPS

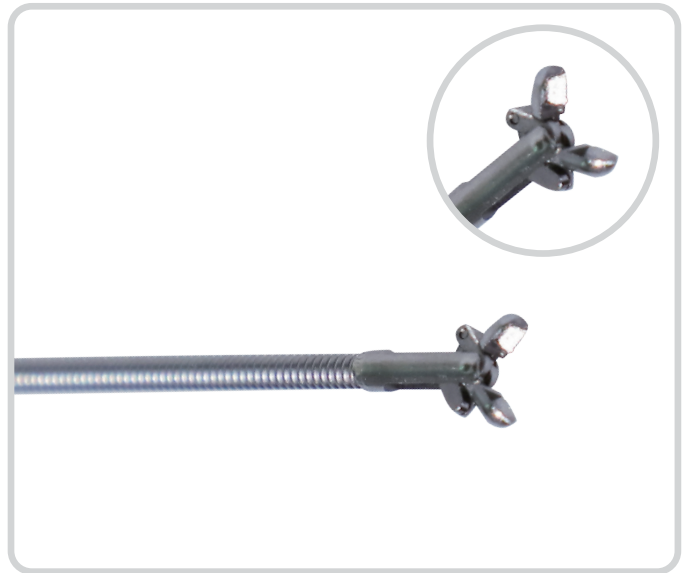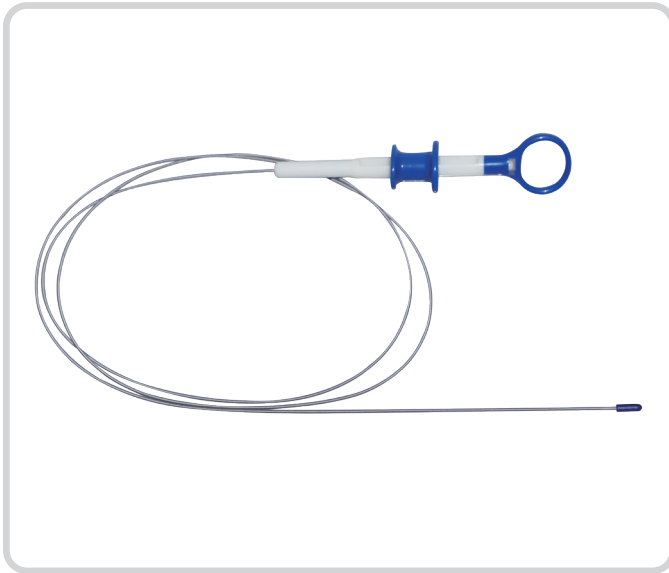

PRODUCT INTRODUCTION

Specification	Features	Purpose
SMH-A 1.8 × 1100	Flat mouth, no needle, no plastic coating on the outer tube, the closing diameter of the tong head is 1.8mm, and the effective length of the product is 1100mm.	Be used for Φ 2.0 and above forceps, clinical use for respiratory tract.
SMH-C+E 1.8 × 1100	Tooth mouth, without needle, plastic coated outer tube, closing diameter of tong head 1.8mm, effective length of product 1100mm.	
SMH-A 2.3 × 1600	Flat mouth, no needle, no plastic coating on the outer tube, the closing diameter of the tong head is 2.3mm, and the effective length of the product is 1600mm.	Be used for Φ 2.8 and above forceps, clinical use for gastroscop.
SMH-B+E 2.3 × 1600	Flat mouth, including needle, plastic coated outer tube, closing diameter of tong head 2.3mm, effective length of product 1600mm.	
SMH-B+E 2.3 × 1800	Tooth mouth, no needle, no plastic coating on the outer tube, the closing diameter of tong head is 2.3mm, and the effective length of the product is 1800mm.	Be used for Φ 2.8 and above forceps, clinical use for 12 fat enteroscopy.
SMH-D+E 2.3 × 1800	Tooth mouth, including needle, plastic coating of outer tube, closing diameter of tong head 2.3mm, effective length of product 1800mm.	
SMH-A 2.3 × 2300	Flat mouth, no needle, no plastic coating on the outer tube, the closing diameter of the tong head is 2.3mm, and the effective length of the product is 2300mm.	Be used for Φ 2.8 and above forceps, clinical use for small enteroscopy.
SMH-D+E 2.3 × 2300	Tooth mouth, including needle, plastic coating of outer tube, closing diameter of tong head 2.3mm, effective length of product 2300mm.	